**gemini**  
ITALY

**BODY COMPOSITION  
MONITORING SCALE**

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MODEL NO 型號: GWS2V1


A GILMAN GROUP BRAND

# IMPORTANT SAFEGUARDS

- **READ ALL INSTRUCTIONS BEFORE USE**
- **PREGNANT WOMAN AND USERS WITH A PACEMAKER OR ELECTRIC EQUIPMENT INSIDE THE BODY SHALL NOT USE THE SCALE.**
- The glass surface of the scale should be kept dry to prevent slipping.
- The scale can only bear normal static pressure when it is weighing. Under abnormal use or impact, there maybe risk of damage to the scale.
- Clean the scale with a damp cloth but do not allow water to go inside the scale and do not use any chemical cleaning agents.
- Do not store the scale in an upright position when it is not in use as this will drain the batteries.
- Remove the batteries if the scale is not being used for a long time.
- The scale is for family use to measure weight, body fat, water mass, muscle mass, bone mass, calorie, BMI. It is not used for commercial or other weighing purpose.
- The body fat parameter is only for your reference. If you target to loss or gain weight, consult your doctor for medical advice.
- In order to get the accurate parameter, please make sure the surrounding area of scale has no electromagnetic objects to avoid interference.
- Please do not repair this scale by yourself. If the scale does not operate normally, check the batteries condition and if the connection direction is right. If the problem cannot be solved, please contact our after-sales service center for assistance.

# PARTS IDENTIFICATION

Before using the scale for the first time, please be familiar with all of the parts. Read all Instructions and Safeguards carefully.

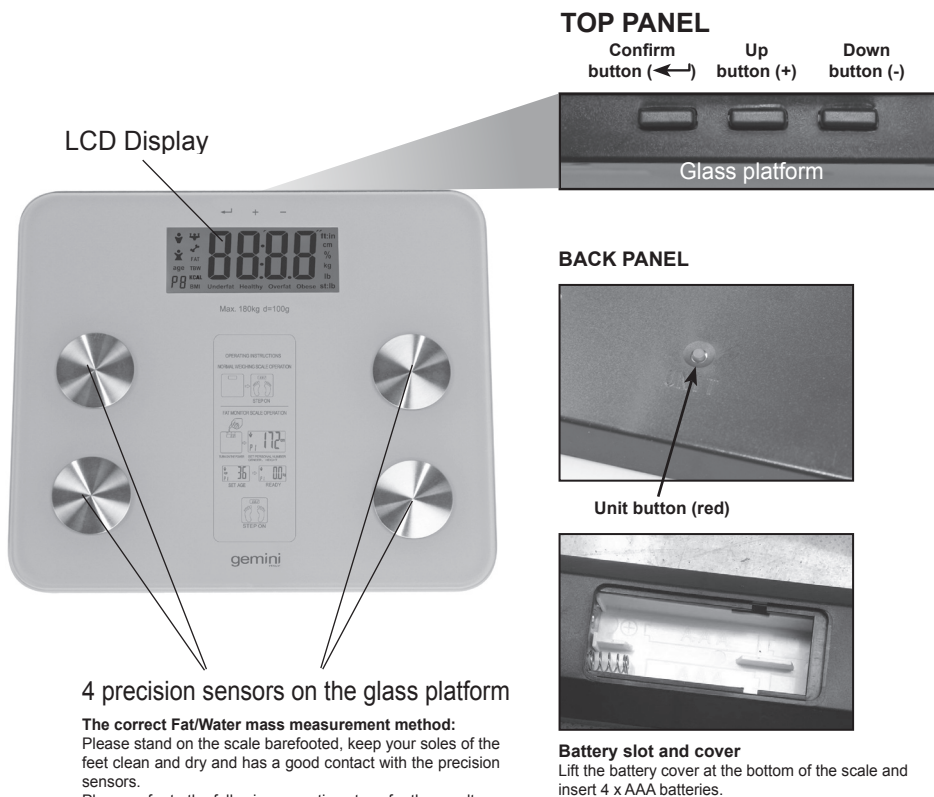4 precision sensors on the glass platform

**The correct Fat/Water mass measurement method:**  
Please stand on the scale barefooted, keep your soles of the feet clean and dry and has a good contact with the precision sensors.  
Please refer to the following operation steps for the result.


## USING THE APPLIANCE

1. Put the scale on a hard and flat floor, avoid to use it on carpet or other soft surface, otherwise it will influence the accuracy.
2. To select the measuring unit at the beginning, press the "Unit" button (red) which is located at the scale's base, the corresponding unit can be selected on the LCD screen, i.e. Kilograms (kg), Stone (st) or Pounds (lb). Choose the desired "Unit" to start the measurement.
3. Stand on the scale directly when using the weight function only.
4. Stand static and put your feet evenly on the glass platform of the scale and wait for the weight calculating result, the correct weight will be displayed on the LCD screen.

**Tips:** The automatic switch on function is only activated for weight function only. You have to input your personal parameters to measure the weight, body fat, water mass, muscle mass, bone mass, calorie and BMI.

# MEASURING WEIGHT, BODY FAT, WATER MASS, MUSCLE MASS, BONE MASS, CALORIE AND BMI

## LCD DISPLAY DETAILS

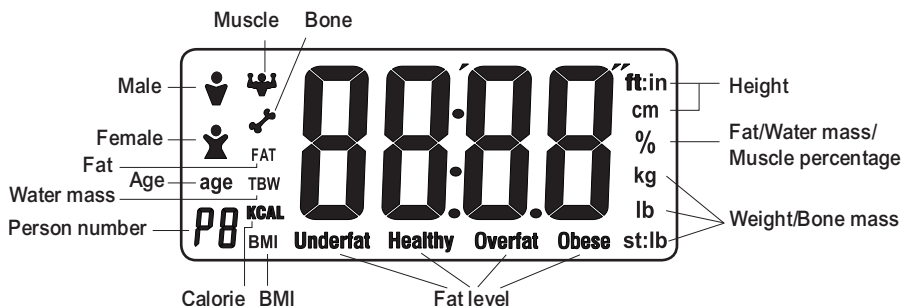1. Put the scale on a hard and flat floor, avoid to use it on carpet or other soft surface, otherwise it will influence the accuracy.
2. Press the "Unit" button (red) at the back of scale, the desired weight unit can be selected, i.e. Kilograms (kg), Stone (st) or Pounds (lb). The corresponding unit will be displayed on the right side of the LCD.
3.
  - a. Press the "Confirm (←→)" button to switch on the screen and wait for "0.0" to appear in display window.
  - b. Select a user number and then input your personal data.

### Enter User Number

- P0 will flash on the LCD display window.
- Press "UP(+)" or "DOWN(-)" button to select your individual user number. (Total 10 users from P0-P9)
- When the selected user number appears on display, press "Confirm (←→)" button to save user number.

### Set Gender

- The gender icon will flash on the LCD display.
- Press "UP(+)" or "DOWN(-)" button to select male (♂) or female (♀) icon.
- Press "Confirm (←→)" button to save gender information.

### Set Height


- The height setting will flash on the LCD display.
- Press “UP(+) or DOWN(-)” button to scroll to your height measurement.
- Press “Confirm (←)” button to save height information.
- The scale is now ready for you to enter your age.

### Set Age

- The “Age” icon will flash on the LCD display.
- Press “UP(+) or DOWN(-)” button to scroll to your age.
- Press “Confirm (←)” button to save age information.
- After the setting, the display will show 0:0 under the preset user number, it is ready to weight.

4. Please make sure your feet are dry and clean when stand on the glass platform, and step on the 4 precision sensors on the glass platform to start weighing, the display will show your weight, body fat, water mass, muscle mass, bone mass, calorie and BMI. These data will be displayed respectively for three times. Each data will display for 2 seconds and then the scale will be switched off automatically.

### Example:

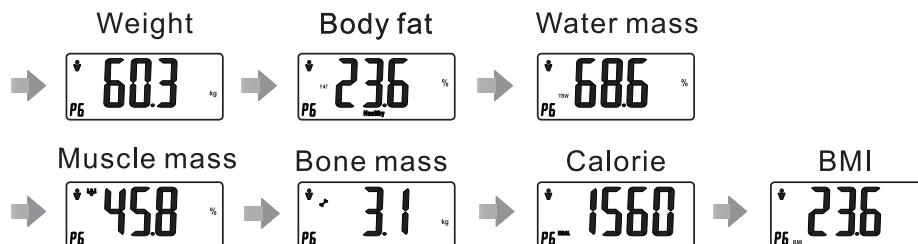5. Please remember your User Number which was set in the beginning. When you use the scale next time, you just need to select your User Number and wait for about 7 seconds, it will display 0.0 on screen. Then stand on the scale to measure your weight, body fat, water mass, muscle mass, bone mass, calorie and BMI.

**Tips:** Please make sure your shoes and socks are taken off before standing on the scale.

**Tips:** When input the data, please do not stand on the glass platform.

## OTHER OPERATION INSTRUCTIONS

1. Fat percentage error indication is **Errr** on screen.
2. Fat percentage lower than 5%, will show **FLo** on screen.
3. Fat percentage higher than 50%, will show **FHI** on screen.
4. Low battery indication is **Lo** on screen.

**\* If you see any of the above error codes, step down from the scale and let the scale off. Try again after 30 seconds.**

# INTERPRETING THE MEASUREMENT DATA

Age	Level of Fat	Male					Female				
		Body fat (%)	Water mass (%)	Muscle mass (%)	Calorie mass	Bone mass (%)	Body fat (%)	Water mass (%)	Muscle mass (%)	Calorie mass	Bone mass (%)
<20	Underfat	<12					<17				
	Healthy	12-20	46-69	38-56	1200-2600	2.0-4.1	17-25	43-68	35-56	1100-2400	1.8-3.9
	Overfat	21-29					26-34				
	Obese	>29					>34				
20-30	Underfat	<13					<17				
	Healthy	13-21	48-69	43-64	1150-2800	2.1-4.0	17-25	43-67	38-58	1050-2600	1.8-3.8
	Overfat	22-30					26-34				
	Obese	>30					>34				
31-40	Underfat	<15					<19				
	Healthy	15-23	45-68	48-69	1100-2400	1.8-4.0	19-27	40-69	42-58	1000-2200	1.5-3.8
	Overfat	24-32					28-36				
	Obese	>32					>36				
41-50	Underfat	<16					<20				
	Healthy	16-24	43-69	46-62	1030-2150	1.9-3.8	20-28	42-68	40-56	950-2050	1.6-3.7
	Overfat	25-33					29-37				
	Obese	>33					>37				
51-60	Underfat	<17					<21				
	Healthy	17-25	43-65	44-59	980-2060	1.9-3.7	21-29	42-65	39-55	950-2000	1.5-3.6
	Overfat	26-34					30-38				
	Obese	>34					>38				
>60	Underfat	<18					<21				
	Healthy	18-26	42-68	41-58	720-1750	1.6-3.7	21-29	41-64	35-51	600-1500	1.3-3.5
	Overfat	27-35					30-38				
	Obese	>35					>38				

The above level of fat is calculated based on the result of body fat (%) and can be seen on screen. The above Water mass / Muscle mass / Bone mass (%) and Calorie show the general range of a Healthy level of fat in terms of Male and Female.

(The above information is for reference only. Please consult your doctor for details and professional advice.)

# BMI STANDARD INDEX FOR DIFFERENT GENDER

1.The BMI (Body Mass Index) is defined as the body weight divided by the square of the body height, and is universally expressed in units of kg/m<sup>2</sup>, resulting from weight in kilograms and height in metres.

$$BMI = \frac{W}{h^2}$$

*w= weight; unit: kilogram*

*h= height; unit: metre*

Gender \ BMI	Under weight	Normal weight	Over weight	Obese
Male	<20	20-25	26-30	>31
Female	<19	19-24	25-29	>30

2. Since BMI does not consider body fat content and other body composition ratios, the BMI range is only a reference based on body weight & height.

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