gemini
ITALY

WEIGHTING SCALE

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MODEL NO 型號: GWS3TS


A GILMAN GROUP BRAND

IMPORTANT SAFEGUARDS

- **READ ALL INSTRUCTIONS BEFORE USE**
- **PREGNANT WOMAN AND USERS WITH A PACEMAKER OR ELECTRIC EQUIPMENT INSIDE THE BODY SHALL NOT USE THE SCALE.**
- The glass surface of the scale should be kept dry to prevent slipping.
- The scale can only bear normal static pressure when it is weighing. Under abnormal use or impact, there maybe risk of damage to the scale.
- Clean the scale with a damp cloth but do not allow water to go inside the scale and do not use any chemical cleaning agents.
- Do not store the scale in an upright position when it is not in use as this will drain the batteries.
- Remove the batteries if the scale is not being used for a long time.
- The scale is for family use to measure weight, body fat, water mass, muscle mass, bone mass, calorie, BMI. It is not used for commercial or other weighing purpose.
- The body fat parameter is only for your reference. If you target to loss or gain weight, consult your doctor for medical advice.
- In order to get the accurate parameter, please make sure the surrounding area of scale has no electromagnetic objects to avoid interference.
- Please do not repair this scale by yourself. If the scale does not operate normally, check the batteries condition and if the connection direction is right. If the problem cannot be solved, please contact our after-sales service center for assistance.

PARTS IDENTIFICATION

Before using the scale for the first time, please be familiar with all of the parts. Read all Instructions and Safeguards carefully.

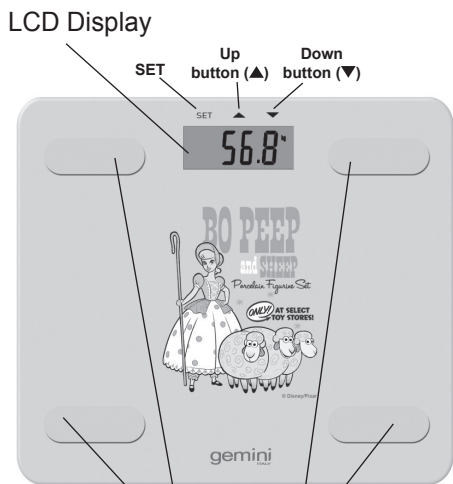
4 precision sensors on the glass platform

The correct Fat/Water mass measurement method:


Please stand on the scale barefooted, keep your soles of the feet clean and dry and has a good contact with the precision sensors.

Please refer to the following operation steps for the result.

BACK PANEL

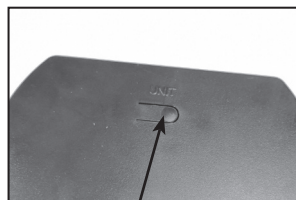Unit button

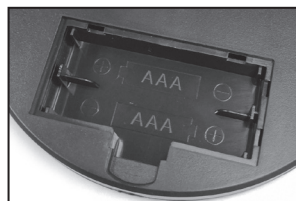Battery slot and cover

Lift the battery cover at the bottom of the scale and insert 2 x AAA batteries.

USING THE APPLIANCE

1. Put the scale on a hard and flat floor, avoid to use it on carpet or other soft surface, otherwise it will influence the accuracy.
2. To select the measuring unit at the beginning, press the “Unit” button which is located at the scale’s base, the corresponding unit can be selected on the LCD screen, i.e. Kilograms (kg), Stone (st) or Pounds (lb). Choose the desired “Unit” to start the measurement.
3. Stand on the scale directly when using the weight function only.
4. Stand static and put your feet evenly on the glass platform of the scale and wait for the weight calculating result, the correct weight will be displayed on the LCD screen.

Tips: The automatic switch on function is only activated for weight function only. You have to input your personal parameters to measure the weight, body fat, water mass, muscle mass, bone mass, calorie and BMI.

MEASURING WEIGHT, BODY FAT, WATER MASS, MUSCLE MASS, BONE MASS, CALORIE AND BMI

LCD DISPLAY DETAILS

1. Put the scale on a hard and flat floor, avoid to use it on carpet or other soft surface, otherwise it will influence the accuracy.
2. Press the “Unit” button at the back of scale, the desired weight unit can be selected, i.e. Kilograms (kg), Stone (st) or Pounds (lb). The corresponding unit will be displayed on the right side of the LCD.
3.
 - a. Press the “SET” button to switch on the screen and wait for “0.0” to appear in display window.
 - b. Select a user number and then input your personal data.

Enter User Number

- P0 will flash on the LCD display window.
- Press “UP(▲)” or “DOWN(▼)” button to select your individual user number. (Total 10 users from P0-P9)
- When the selected user number appears on display, press “SET” button to save user number.

Set Gender

- The gender icon will flash on the LCD display.
- Press “UP(▲)” or “DOWN(▼)” button to select male (♂) or female (♀) icon. For some users who exercises regularly, please select Sport male (♂) (♂) or Sport female (♀) (♀).
- Press “SET” button to save gender information.

Set Height

- The height setting will flash on the LCD display.
- Press “UP(▲)” or “DOWN(▼)” button to scroll to your height measurement.
- Press “SET” button to save height information.
- The scale is now ready for you to enter your age.


Note: When the measurement unit is set to kilogram (kg), the unit of height is centimeter (cm); When the measurement unit is set to stone (st) or pound (lb), the unit of height is inches (’ ’’).

Set Age

- The “Age” icon will flash on the LCD display.
- Press “UP(▲)” or “DOWN(▼)” button to scroll to your age.
- Press “SET” button to save age information.
- After the setting, the display will show 0:0 under the preset user number, it is ready to weight.

4. Please make sure your feet are dry and clean when stand on the glass platform and step on 4 precision sensors on the glass platform to start weighing, the 4“0” display one by one from left to right, then the display will show your body fat, water mass, BMI, muscle mass, bone mass, calorie and weight. These data will be displayed respectively for two times.

Example:

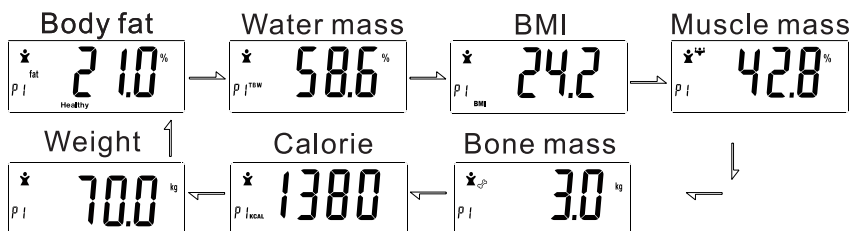5. Please remember your User Number which was set in the beginning. When you use the scale next time, you just need to select your User Number and wait for about 3 seconds, it will display 0.0 on screen. Then stand on the scale to measure your weight, body fat, water mass, muscle mass, bone mass, calorie and BMI.

Tips: Please make sure your shoes and socks are taken off before standing on the scale.

Tips: When input the data, please do not stand on the glass platform.

OTHER OPERATION INSTRUCTIONS

1. Overload indication is **Err** on screen.
2. Low battery indication is **Lo** on screen.

*** If you see any of the above error codes, step down from the scale and let the scale off. Try again after 10 seconds.**

INTERPRETING THE MEASUREMENT DATA

Body fat analysis (Unit %)								
Age	Female				Male			
	Underfat	Healthy	Overfat	Obese	Underfat	Healthy	Overfat	Obese
10-39	<21	21-33	33-39	>39	<8	8-20	20-25	>25
40-59	<23	23-34	34-40	>40	<11	11-22	22-28	>28
60-80	<24	24-36	36-42	>42	<13	13-25	25-30	>30

Body water analysis (Unit %)						
Age	Female			Male		
	Underfat	Healthy	Overfat	Underfat	Healthy	Overfat
10-15	<57	57-67	>67	<58	58-72	>72
16-30	<47	47-57	>57	<53	53-67	>67
31-60	<42	42-52	>52	<47	47-61	>61
61-80	<37	37-47	>47	<42	42-56	>56

Body muscle analysis (Unit %)		
Gender	Female	Male
Normal	>34%	>40%

Body bone analysis (Unit Kg)						
Gender	Female			Male		
Weight	<45Kg	45-60Kg	>60Kg	<60Kg	60-75Kg	>75Kg
Bone mass	1.8Kg	2.2Kg	2.5Kg	2.5Kg	2.9Kg	3.2Kg

The above level of fat is calculated based on the result of body fat (%) and can be seen on screen. The above Water mass / Muscle mass / Bone mass (%) and Calorie show the general range of a Healthy level of fat in terms of Male and Female.

(The above information is for reference only. Please consult your doctor for details and professional advice.)

BMI STANDARD INDEX FOR DIFFERENT GENDER

1.The BMI (Body Mass Index) is defined as the body weight divided by the square of the body height, and is universally expressed in units of kg/m², resulting from weight in kilograms and height in metres.

$$BMI = \frac{W}{h^2}$$

w= weight; unit: kilogram

h= height; unit: metre

Gender \ BMI	Under weight	Normal weight	Over weight	Obese
Male	<20	20-25	26-30	>31
Female	<19	19-24	25-29	>30

2. Since BMI does not consider body fat content and other body composition ratios, the BMI range is only a reference based on body weight & height.

TECHNICAL DATA

Item no: GWS3TS

Power: 2 x AAA batteries (Included)

Product Dimension: 300(W) x 22(H) x 265(D)mm

Maximum Capacity: 396 lb / 180 kg / 28.34 ST

* The final product specification may be varied according to a different lot of material and actual input voltage power of the usage area.

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